



TXD3 CYTO Centrifuge

TXD3 cell smear centrifuge can be used for the liquid handling in fields of cytology, hematology and microbiology. Such as cerebrospinal fluid, pleural fluid, peritoneal fluid, synovial fluid, urine, oral cavity-cleaning, stomach cleaning, exudate, percolate, phlegm. It is mainly used for checks with small cell capacity. Our clip rotor for cytology is compatible with popular cytocentrifuge accessories currently on the market.

Product advantages

- Microprocessor control and DC brushless motor.
- Touch panel and LCD display. RCF can be calculated automatically.
- Save up to 25 user-defined programs. Featuring a separate short spin key for quick centrifugation.
- Special damping structure to reduce vibration.
- Centrifuge lid with soft-touch lid closure. centrifuge cannot work if the door is open and the door cannot be opened when it is working.
- Several smear rotors available, it can smear various fluids to achieve the best centrifugation effect.

Technical parameter

Order NO	0205198000
Max Speed	3000r/min
Max RCF	890xg
Funnel Capacity	0.1ml-0.5ml
Rotor Capacity	12pcs
Speed Accuracy	±30r/min
Noise	≤ 60dB (A)
Power Supply	AC 220±22V 50/60Hz
Power Consumption	260W
Dimension (WxDxH)	380×500×300(mm)
Packing Dimension (WxDxH)	465×660×350 (mm)
Net Weight	19kg

Rotor parameters



Swing bucket cell smear rotor
Max Capacity: 6pcs
Max Speed: 2000r/min
Max RCF: 545xg
Order NO: 0302855810



Bucket
Order NO: 0303855804



Plastic holder
Order NO: 0304855808



Cell smear rotor
Max Capacity: 6pcs
Max Speed: 3000r/min
Max RCF: 730xg
Order NO: 0302856510



Smear clip
Order NO: 0308855521




Glass slide
Order NO: 0310856514



Cell smear rotor
Max Capacity: 12pcs
Max Speed: 3000r/min
Max RCF: 890xg
Order NO: 0302856520



Filter paper
Order NO: 0308855524

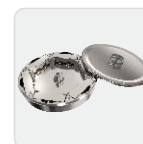

Bio-seal smear rotor
Max Capacity: 6 Plastic holder
Max Speed: 2000r/min
Max RCF: 405xg
Order NO: 0301855530